

The dynamically crash-tested **PRO F1 Universal Cot Fastener** is the first truly SAE-compliant cross-compatible fastening system. It is designed for use with:







Stryker Power-PRO XT



Stryker Performance-PRO XT

With the PRO F1, you can easily exchange cots between the vehicles in your fleet as needed, knowing that they will be safely secured in compliance with crash-test standards.



One-touch foot-end release frees cot from the PRO F1 with ease.



SAE J3027-compliant for patient and provider safety in the ambulance.



Cross-compatible design allows for use with FERNO and Stryker cot models.






- Universal design accommodates FERNO POWER X1[™], Stryker Power-PRO[™] XT, and Stryker® Performance-PRO™XT.
- Securely fastens compatible cots during transport for patient and provider safety
- · Linear guide system facilitates loading and unloading, even on an angle, reducing lateral cot drift
- · Reduced lateral movement allows for more efficient loading and unloading



Specifications

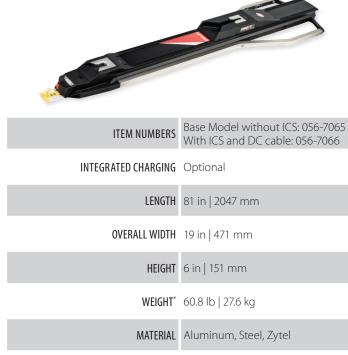

COT

Compatible Systems

STANDARD

Features & Benefits

- Designed for longevity with durable aluminum, steel, and Zytel® construction
- Low profile and enclosed construction minimizes trip hazards
- One-touch foot-end release button easily releases cot from fastener
- Easy to remove and power washable for cleaning and decontamination
- 110° approach angle when combined with the POWER X1 Ambulance Cot
- ICS® Integrated Charging System integration available for use with the FERNO POWER X1 Ambulance Cot (features a remote charging indicator)
- · Compatible with the FERNO Universal Floor Plate
- 2-Year Warranty

Compliance

The **PRO F1** has been designed to meet international test standards to increase safety for you and your patient.

- SAE J3027
- KKK-A-1822F
- IEC 60601-1-2

- CAAS-GVS v1.0
- NFPA 1917
- ISTA 1A

• IP56

COT FASTENER

- AMD 004
- EN 1789

DATIENT DECTRAINT

The PRO F1 is compliant with the above standards when used with:

- POWER X1 combined with the PRO F1 Compatibility Kit
- Stryker Power-PRO XT configured for the Stryker Performance-LOAD or Stryker Power-LOAD***
- Stryker Performance-PRO XT configured for the Stryker Performance-LOAD or Stryker Power-LOAD***

STANDAKD	=	COI		+ COLFASIENER				+ PATIENT RESTRAINT			
	FERNO POWER X1	Stryker Performance-PRO XT	Stryker Power-PRO XT	FERNO PRO F1	Stryker Performance-LOAD	FERNO Rail/Antler Series	Stryker Rail/Antler Series**	FERNO X-Restraints	FERNO BioSafe Restraints	Stryker X-Restraints	Stryker Rugged-X Restraints**
	•	Teriorinance Prio XI	TowerThoxi	•	Terrormance cond	null/Alltici Sches	nui/Airaci Scries	• Nestraines	biosaic nestraints	x nestraines	nestraines
	•			•					•		
	•				•			•			
		•		•						•	
			•	•						•	
		•		•							•
			•	•							•
CAAS-GVS v1.0 (2016) (SAE J3027)	•			•				•			
	•			•				•	•		
		•			•			•			
		•	•								
		•		•							•
			•	•							•
NFPA 1917 (2016) (SAE J3027)	•			•				•			
	•			•					•		
	•				•			•			
		•		•						•	
			•	•						•	
		•		•							•
			•	•							•
KKK-A-1822F through change notice 8 (AMD 004)	•					•		•			
	•			-		•		•	•		
							•	•	•		
						•	•		-		
		<u> </u>									+

^{**} As advertised

When used in conjunction with the manufacturer supplied floor plate system and safety hook, each of the combinations listed in the matrix complies with the respective standard listed.

""Power-PRO XT, Performance-PRO XT, Performance-LOAD*, and Power-LOAD* are trademarked products of Stryker Corporation. FERNO is not associated with Stryker Corporation. FERNO reserves the right to change specifications without notice.

^{*} Fastener is weighed without ICS and mounting hardware